



Ultra-Viol sp.j. Pietras, Purgał, Wójcik Stępowizna 34, 95-100 Zgierz, Poland VAT PL 727 00 21 903 www.ultraviol.pl/en Contact: phone number +48 42 717 77 45 ext. 21 mobile +48 725 779 277 fax +48 42 715 02 16, e-mail: office@ultraviol.pl


TECHNICAL DATA SHEET **NBV series** 



### **NBV 15 N** DIRECT RADIATION UV-C GERMICIDAL LAMP

# HOW DOES THE DIRECT RADIATION UV-C GERMICIDAL LAMP WORK?

Direct radiation germicidal lamps type NBV are designed to prevent primary and secondary infections of patients and medical personnel caused by airborne pathogenic microorganisms (pathogens). Using direct germicidal radiation in the rooms where infected patients or patients with immune deficiencies are staying, significantly reduces the probability of infection spread by air. Raising the level of microbiological purity of the air and the rooms helps to destroy and reduce impact of existing outbreaks of pathogens.

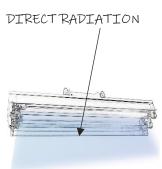

CLEAN AIR AND SURFACE

#### **BASIC DATA:**

no
yes
not allowed
wall-mounted
no
yes, 1 bulb
powder laquered aluminium
high quality reflective aluminium
no
+10°C to +40°C
30% to 70%
700 hPa to 1060 hPa
yes





Ultra-Viol sp.j. Pietras, Purgał, Wójcik Stępowizna 34, 95-100 Zgierz, Poland VAT PL 727 00 21 903 www.ultraviol.pl/en Contact: phone number +48 42 717 77 45 ext. 21 mobile +48 725 779 277 fax +48 42 715 02 16, e-mail: office@ultraviol.pl

### **TECHNICAL DATA SHEET**

# **NBV** series

## NBV 15 N

### DIRECT RADIATION UV-C GERMICIDAL LAMP

### **TECHNICAL DATA:**

Supply voltage	230 V, 50 Hz
Power consumption	18 W
UV-C bulbs (Philips/Osram)	1 x 15 W (TUV/HNS)
UV-C radiation wave-length	253.7 nm
Useful lifetime of the UV-C bulbs	min. 8000 h
Radiation intensity of the external UV-C tube at the distance of 1 m $$	0.9 W/m <sup>2</sup>
Effective area of the lamp	6-8 m <sup>2</sup>
Exposure angle adjustment range	160°
Class of protection against electric shock	I
Ingress Protection Code	IP 20
Operation mode	continuous
Lamp body dimensions - without holder $(L \times W \times H)$	500 x 85 x 135 mm
Total lamp mass	2.4 kg
Holder length	12.0 cm
Power cord length	0.8 m

### **APPLICATION RECOMMENDATIONS:**

#### Hospitals

- operation theatres
- intensive care units
- emergency rooms
- examination and treatment rooms
- reception units
- patient rooms, isolation rooms
- soiled/dirty utility rooms

Outpatient clinics Medical laboratories Chemist's Beauty salons Pharmaceutical industry Food industry Cosmetic industry



Ultra-Viol's realisation