

TECHNICAL DATA SHEET

NBV series




Ultra-Viol sp.j. Pietras, Purgał, Wójcik Stępowizna 34, 95-100 Zgierz, Poland VAT PL 727 00 21 903 www.ultraviol.pl/en Contact: phone number +48 42 717 77 45 ext. 21 mobile +48 725 779 277 fax +48 42 715 02 16, e-mail: office@ultraviol.pl



NBV 2x30 SLW DIRECT RADIATION UV-C GERMICIDAL LAMP

HOW DOES THE DIRECT RADIATION UV-C GERMICIDAL LAMP WORK?

Direct radiation germicidal lamps type NBV are designed to prevent primary and secondary infections of patients and medical personnel caused by airborne pathogenic microorganisms (pathogens). Using direct germicidal radiation in the rooms where infected patients or patients with immune deficiencies are staying, significantly reduces the probability of infection spread by air. Raising the level of microbiological purity of the air and the rooms helps to destroy and reduce impact of existing outbreaks of pathogens.

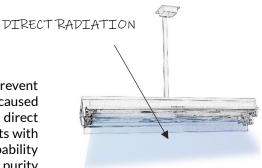

CLEAN AIR AND SURFACE

BASIC DATA:

Air flow disinfection function (flow disinfection chamber)	no
Direct radiation function	yes
Presence of people, animals, plants inside the room during disinfection process	not allowed
Mounting type	ceiling-mounted
Working time counter	yes, digital working time counter with 4-field LED display and acoustic signalling
External bulbs	yes, 2 bulbs
Casing material	powder laquered aluminium
Reflector material	high quality reflective aluminium
Switch on/off	no
Ambient temperature	+10°C to +40°C
Relative humidity	30% to 70%
Atmospheric pressure	700 hPa to 1060 hPa
Declaration of conformity	yes
User's manual English version	yes





Ultra-Viol sp.j. Pietras, Purgał, Wójcik Stępowizna 34, 95-100 Zgierz, Poland VAT PL 727 00 21 903 www.ultraviol.pl/en Contact: phone number +48 42 717 77 45 ext. 21 mobile +48 725 779 277 fax +48 42 715 02 16, e-mail: office@ultraviol.pl

TECHNICAL DATA SHEET

NBV series

NBV 2x30 SLW

DIRECT RADIATION UV-C GERMICIDAL LAMP

TECHNICAL DATA:

Supply voltage	230 V, 50 Hz
Power consumption	66 W
UV-C bulbs (Philips/Osram)	2 x 30 W (TUV/HNS)
UV-C radiation wave-length	253.7 nm
Useful lifetime of the UV-C bulbs	min. 8000 h
Radiation intensity of the external UV-C tube at the distance of $1\mathrm{m}$	3.6 W/m ²
Effective area of the lamp	18-22 m ²
Exposure angle adjustment range	240°
Class of protection against electric shock	I
Ingress Protection Code	IP 20
Lamp body dimensions - without holder $(L \times W \times H)$	960 x 85 x 145 mm
Total lamp mass	4.6 kg
Extension arm (holder) length:	50.0 cm
Power cord length	0.8 m

APPLICATION RECOMMENDATIONS:

Hospitals

- operation theatres
- intensive care units
- emergency rooms
- examination and treatment rooms
- reception units
- patient rooms, isolation rooms
- soiled/dirty utility rooms

Outpatient clinics Medical laboratories Chemist's Beauty salons Pharmaceutical industry Food industry Cosmetic industry



Ultra-Viol's realisation